

Using GC-MS/MS for Superior Sensitivity, Specificity and Precision in Free Testosterone Analysis

This article discusses biologically important steroid measurements which have, in the past, seen a lack of analytical sensitivity. Examples include determination of testosterone levels in women and children and oestrogen levels in post-menopausal women. The article gives insight into the analysis of free testosterone at physiologically relevant levels in plasma ultra-filtrate, concentrating on the sensitivity, specificity, precision and accuracy advantages offered by GC-MS/MS.

Although highly efficient for the analysis of complex mixtures, GC-MS/MS could not address the requirement of the industry for the analysis of larger non-volatile molecules...

Author Details:

By Rohan A. Thakur, Clark Williard and Agnita Rajasekaran
Taylor Technology - a subsidiary of PharmaNet Development Group, Princeton, NJ

INTRODUCTION

Until 1955, mass spectrometry (MS) was most commonly used for the direct analysis of volatiles. Gas chromatography (GC) was coupled to MS for the first time at Dow Chemical Company with the aim to expand the analytical capabilities of MS to cover complex mixtures of unknowns [1]. Subsequent technological developments involved the introduction of hybrid mass analyzers such as the triple stage quadrupole (TSQ) mass spectrometers and the use of tandem MS (MS/MS or SRM) as a high specificity technique for routine quantitative analysis of complex mixtures. GC-MS/MS was quickly established as the technology of choice for bioanalytical applications.

Although highly efficient for the analysis of complex mixtures, GC-MS/MS could not address the requirement of the industry for the analysis of larger non-volatile molecules. To overcome this new challenge, high performance liquid chromatography (HPLC) was coupled to MS/MS (LC-MS/MS) leading to the development of the atmospheric pressure based ionisation (API) technique of electrospray (ESI). This new method soon was adopted as the preferred technology for bioanalysis [2].

In addition, LC-MS/MS was seen as offering a number of advantages that allowed for easy, uncomplicated bioanalysis of small molecules. Compared to GC-MS/MS, the technique requires very little sample preparation while delivering good detection limits and is also capable of analysing a wider range of compound classes and molecular weights. Latest advances have seen the introduction of TurboFlow, on-line solid phase extraction (SPE) and other chromatographic techniques that facilitate direct injection of biological fluids into LC-MS/MS systems, eliminating sample preparation [3]. These methods are particularly beneficial for high-throughput laboratories, where rapid analysis of thousands of samples is needed. However, in applications where assay precision and accuracy are of utmost importance, LC-MS/MS may not always deliver the required results. Indeed for the highly sensitive and specific detection of steroids for example, it has been suggested that GC-MS/MS might actually be the preferred technique.

Oestrogens (oestradiol and oestrone) and androgens (testosterone and androstenedione) are a class of biologically important endogenous compounds that routinely require reproducible results at very low detection limits (attogram level). Endogenous and exogenous oestrogens are associated with the pathogenesis of breast cancer [4] whereas testosterone affects major biological functions and is used as a therapeutic agent for the treatment of hypogonadism, the stimulation of erythropoiesis and hereditary angioneurotic edema [5]. As a consequence, measurements of oestrogens and androgens is of high biological importance, determining for example testosterone levels in women and children and estrogen levels in post-menopausal women [6,7]. Analytical sensitivity and specificity play major roles in the interpretation of steroid measurements. Sample preparation is critical for the sensitivity, precision and accuracy of results. Ion suppression can be directly related to inadequate sample preparation and is a major problem with ESI based LC-MS/MS techniques [8].

The limitations of LC-MS/MS related to accuracy, precision, specificity and sensitivity have been significant enough to initiate an inter-laboratory comparison study of serum total testosterone measurements performed by mass spectrometric methods [6]. Recent data indicates that GC-MS/MS offers three key advantages over LC-MS/MS when it comes to this type of analysis. The derivatisation technique improves sensitivity (Neg Ion CI dramatically

reduces chemical noise), the high number of theoretical plates of a GC column results in highly focused chromatographic peaks (approximately 3 seconds or less), and the entire column eluent enters the ion source (unlike LC-MS/MS ion sources which sample 5% of the LC eluent). As a result, GC-MS/MS provides the specificity, sensitivity and precision required to reliably interpret physiologically relevant levels of free testosterone in plasma.

An experiment was developed to demonstrate the high sensitivity, precision and accuracy offered by GC-MS/MS. A bioanalytical method was used for the analysis of free testosterone at physiologically relevant levels in plasma ultra-filtrate (UF) from post-menopausal women.

EXPERIMENTAL

EDTA (ethylenediaminetetraacetic acid) plasma from post-menopausal women was stored at -20°C and used to form quality control samples (QCs). The plasma was also stripped with charcoal to prepare calibration standards (STDs). Plasma UF was prepared using a centrifugal filter device (Centricon YM-30, Millipore Corporation, Billerica, MA). Sample preparation was completed in three steps. First, a penta deuterated testosterone (d5- testosterone) solution was added to 0.5mL of the plasma UF, followed by liquid:liquid hexane/ethyl acetate extraction. Oxime of the steroids was then prepared by evaporating a volatile extraction solvent with N₂ and derivatising with O-(2,3,4,5,6 pentafluorobenzyl)-hydroxylamine HCl. Silylation of the OH group was enabled using N-Methyl-N-trifluoroacetamide (MSTFA).

A gas chromatograph (Varian 3400, Varian, Inc., Walnut Creek, CA) was implemented together with a 15m x 0.250mm DB-1 fused silica capillary column (J&W Scientific, Folsom, CA). The split flow rate was 50:1 and the temperature program spanned from 210°C to 303°C. A triple stage mass spectrometer (Finnigan TSQ-7000, Thermo Fisher Scientific, San Jose, CA) was also used, operated in negative ion chemical ionisation (NICI) mode.

RESULTS AND DISCUSSION

The assay range of the GC-MS/MS method was determined to be between 0.250 and approximately 100 pg/mL. The calibration standard curves yielded consistent regression slopes and all acceptance criteria were met. *Table 1* demonstrates the intrarun and interrun precision and accuracy results whereas *Figure 1* shows the chromatograms that were derived over the course of the experiment. Recovery was 78-79%. UF samples were stable after three freeze/thaw cycles, after 30.5 hours at ambient temperature, after 42 days at -20°C and after 46 days at -70°C. Extracts were stable in the autosampler after 239 hours. Testosterone levels were stable in human plasma stored at -70°C for at least 62 days and -20°C for at least 42 days.

GC-MS/MS was able to routinely and reproducibly deliver results at a LOQ of 3 femtograms on column, considerably lower than most LC-MS/MS systems can routinely deliver. The liquid-liquid extraction technique used for this assay is proven to reduce matrix effects, resulting in the analytical column always being exposed to the cleanest possible biological samples. As a result, build-up of problematic matrix components such as phospholipids, polar metabolites and other polar interferences was reduced. (8) In addition, the derivatisation step required to make compounds volatile for injection into the GC offered some degree of specificity.

The GC column provided increased efficiency of chromatographic separation, generating highly focused peaks. Since the majority of the analyte was chromatographically concentrated in a narrow peak (3

seconds or less at base), separation from background noise became easier, resulting in clear delineation as to when the signal started and ended. Consequently, precision was improved since peak inflection and deflection points were easier to assess through automated peak integration algorithms.

The entire column effluent from the GC column entered the MS ion source making it far more efficient than atmospheric pressure ionisation (API) techniques employed by LC-MS/MS. API sources (ESI/APCI) are usually 1-10% efficient in transferring the ions from the LC eluent into the MS. This means while close to 100% of the column eluent enters the GC-MS/MS ion source, exactly the opposite occurs in an LC-MS/MS API ion source. Until this sampling rate is improved, it appears that GC-MS/MS will be more efficient and this has a direct impact on sensitivity.

CONCLUSION

Choosing whether to implement LC-MS/MS or GC-MS/MS for bioanalysis depends on the specific application requirements. LC-MS/MS is the preferred technique for simple bioanalysis of an extended range of compound classes and molecular weights from small to large non-volatile molecules. However, the method does not offer the accuracy, precision, specificity and sensitivity required for specific applications including the analysis of steroids (estrogens and androgens). A highly sensitive GC-MS/MS bioanalytical assay for free testosterone in UF plasma from post-menopausal women has been successfully implemented to help elucidate the pharmacokinetic parameters of testosterone-related therapy and the monitoring of endogenous testosterone as a pharmacodynamic biomarker.


REFERENCES

- Gohlke, R.S. and McLafferty, F. W. J. *Amer Soc. Mass Spectrom.* (1993), 4, 367
- Yamashita, M, and Fenn, J.B. *J. Phys. Chem.* (1984) 88, 4671.
- Herman J. L. *Rapid Commun. Mass Spectrom.* (2002) 16, 421.
- Wang, S, Paris, F, Sultan, C. S, Song, R., Demers, L.N, Sudaram, B., Settlege, J, Ohorodnik, S, and Santen, R. *JCEM* (2005) 3, 1407.
- Santen, R, Demers, L, Ohorodnik, S, Settlege, J, Langecker, P, Blanchett, D, Goss, P. E, and Wang S. *Steroids* (2007) 72, 656.
- Vesper, H.W, Bhasin, S, Wang, C, Tai, S, Dodge, L.A, Singh, R. J, Nelson, J., Ohorodnik, S, Clarke, N. J, Salameh, W. A, Parker, R. C. Jr, Razdan, R., Monsell, E. A. and Myers, G. L. *Steroids* (2009) 74, 498.
- Singh. R. J. *AACC* (2008) 34, 12.
- Bonfiglio, R, King, R.C, Olah, T. V, Merkle, K. *Rapid Commun. Mass Spectrom.* (1999) 13, 1175.

We would like to acknowledge Robert White, Randy Weintraub, Bhaskar Sundaram, James Settlege, and Susan Ohorodnik.

Table 1. Intra-run/Inter-run Precision and Accuracy. LLOQ is the lower limit of quantitation and ULOQ is the upper limit of quantitation.

		LLOQ	Low	Medium	ULOQ
		0.250 pg/mL	1.29 pg/mL	9.56 pg/mL	96.1 pg/mL
Intra-run					
Run 1	Mean	0.242	1.37	9.36	94.9
	CV (%)	14.8	6.69	7.18	4.58
	Bias (%)	-3.2	6.20	-2.09	-1.25
	N	5	6	6	6
Run 2	Mean	0.222	1.25	9.61	96.1
	CV (%)	10.4	3.26	2.26	6.50
	Bias (%)	-11.2	-3.10	0.523	0.00
	N	6	6	6	6
Run 3	Mean	0.207	1.25	9.73	97.3
	CV (%)	22.5	7.00	0.741	2.70
	Bias (%)	-17.2	-3.10	1.78	1.25
	N	5	5	6	6
Inter-run					
	Mean	0.224	1.29	9.56	96.1
	CV (%)	16.1	7.16	4.36	4.66
	Bias (%)	-10.4	0.00	0.00	0.00
	N	16	17	18	18

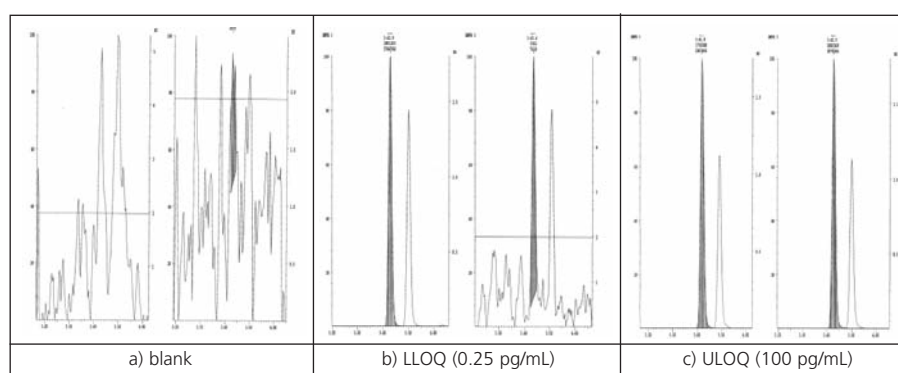

Note: in each panel, the first chromatogram is a trace of the IS ion channel; the second is a trace of the ion channel for testosterone. The second peak in the LLOQ and ULOQ chromatograms is due to the testosterone isomer formed during derivatisation.

Figure 1. Chromatograms of Extracted Standards that are: a) blank, b) at the LLOQ, and c) at the ULOQ.

High Pressure Chemical Reactors

Supercritical Fluid Technologies introduces a new series of stirred reactors for high-pressure chemistry. The HPR Series Reactors have been designed for researchers who are interested in performing pressurised chemical reactions in their laboratories. Whether it is high-pressure chemical synthesis or process development, the HPR Series is an excellent unit for those applications. The reactors range in size from 50ml to 4 litres and may be operated up to 10,000 psi and 350°C. The reactors have a magnetically coupled impeller for optimal mixing. All high-pressure components are ASME rated and protected by a rupture disc for safe operation. These laboratory bench top models are ideal for applications where repetitive use makes convenience a necessity, such as catalytic studies, polymerisation, hydrogenation, oxidation, isomerisation, and dehydrogenation. All sizes are supplied as ready-to-use instruments requiring only utility connections prior to operation. The reactors are compact instruments that fit easily into a fume hood. Their modular design makes it simple and cost-effective to alter the unit's basic configuration to adapt to new or evolving application needs.

All functions for the HPR Series Reactor are controlled by choices of two integrated processors (RxTrol Jr. or RxTrol Sr.) with a full colour touch screen. RxTrol Jr. consists of (2) temperature sensors (wall/interior mounting), and fuzzy logic PID temperature controls and mixing controls in a table top assembly. RxTrol Sr. offers many more advanced features utilising a microcontroller to perform all PID, ramp/soak, process trending, and interlock functions. Additional options include vessel windows (<150°C operation only), additional ports, cooling coils, sample loops, stand alone pressure transducers, supercritical fluids pumps and reagent addition modules.

Circle no. 310



New FTIR Spectral Libraries Aid Identification of Rocks and Minerals



A2 Technologies has announced the availability of two new diffuse reflectance mid-infrared spectral libraries of rocks and minerals further extending the utility of their innovative Exoscan and FlexScan handheld FTIR analysers for geoscience applications. The first library contains spectra of general rocks and minerals sourced from around the world, while the other is focused on commercially important or strategic minerals.

FlexScan and Exoscan are the world's first FTIR analysers that have the performance and robustness to be applicable to field based geoscience investigations. The analysers are equipped with diffuse reflectance sampling interfaces, thus enabling rocks, minerals and soils to be analysed with minimal, or no sample preparation. Now, minerals of interest located in outcroppings, ore bodies, drill cores, and mines can be analysed directly at the locale. The capability of FTIR to elucidate how atoms are bound together to form the chemistry of rocks and minerals is an important complement to widely used hand-held XRF analysers, which identify the elements that comprise the geological sample.

The diffuse reflectance infrared spectra of rocks and minerals are different than those that are contained in classical transmission IR libraries that already exist. For this reason, A2 Technologies is building, and will continue to build, these unique diffuse reflectance libraries in support of its efforts to bring FTIR analysers to geoscience applications. Since the chemistry of rocks and minerals is complex and highly dependent on the specific locale and condition of formation, A2 Technologies analysers permit users to easily create their own diffuse reflectance IR libraries to meet the requirements of their particular area of interest.

Circle no. 311